



*****NEWBORN SCREENING COLLECTION INFORMATION*****

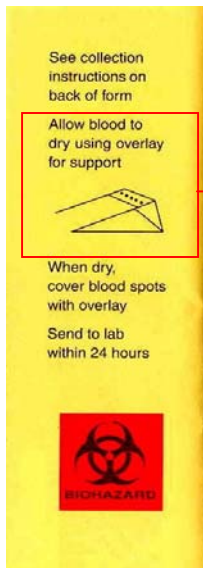


Do not use this drying technique displayed on the overlay flap; too much of a drying angle may cause serum to separate. Blood may also come in contact with edge of flap and be drawn out of filter.

It is very important that the wet filter paper **does not** come in contact with any surface until completely dry. The filter circles are designed to hold a specific quantity of blood, if the wet filter is allowed to come in contact with the gold paper backing of form or any other surface, blood can be drawn out of the filter making the quantitative tests performed by the Newborn Screening Laboratory invalid. The gold colored backing and flap portion can be folded away from the filter for collection and drying of blood specimen.

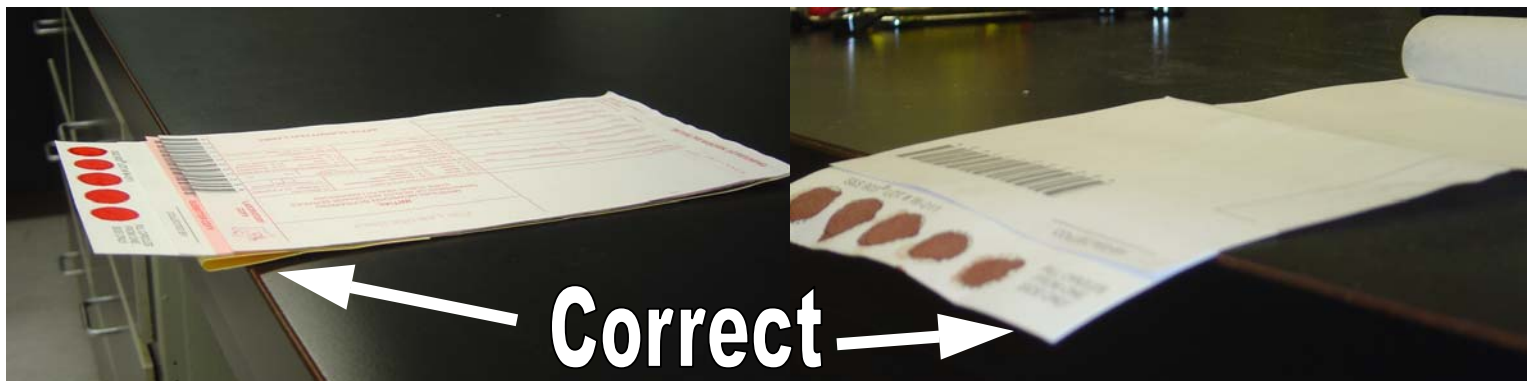

Please do not allow the wet filter paper to rest on the gold colored portion of the form

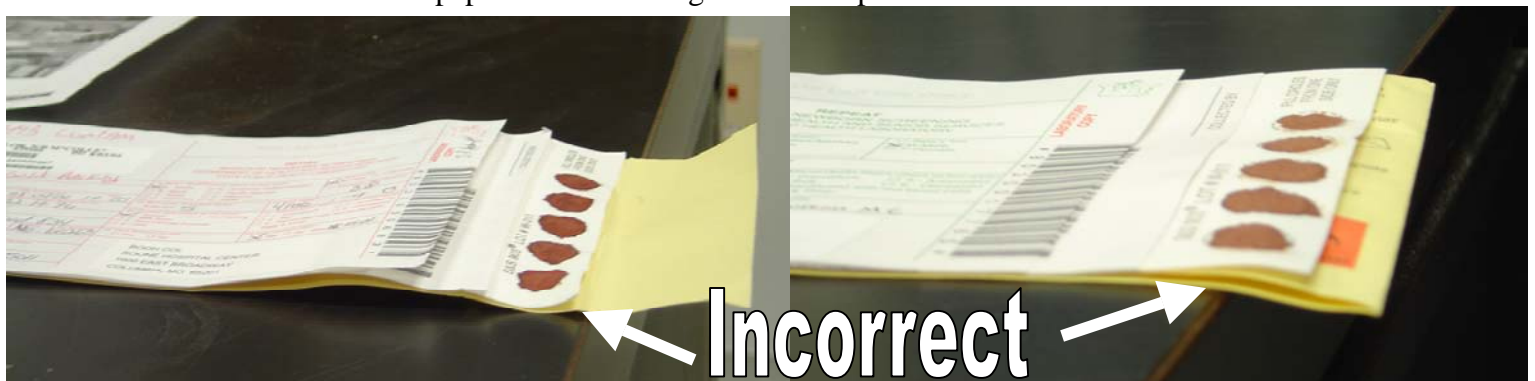

Or blood may be drawn out of the filter and onto the gold colored backing of form.

